





INNOVATIONS IN DENTISTRY SINCE 2000

amdlasers.com

1.866.999.2635

support@amdlasers.com



Picasso Lite Plus

PLPLUS

The Picasso Lite+ Diode Laser was designed for versatility, precision, and affordability. Offering a wide range of capabilities for soft tissue procedures and periodontal treatments—it's the ideal laser for both dentists and hygienists alike.

Picasso Plus

PPLUS

As AMD's top-tier laser, the powerful Picasso+ Diode Laser offers a complete solution for dental professionals looking for power, control, and flexibility. The Plus is the only Picasso that comes with an attenuator tip for laser-assisted teeth whitening.

Picasso Clario

PCPLUS

The Picasso Clario Diode Laser was designed for the hygiene role, complete with streamlined navigation and pre-set parameters based on diagnosis, which facilitates simple, controlled, and easy use of the laser for hygiene treatments and soft tissue procedures.

PICASSO DISPOSABLE TIPS AND ACCESSORIES

Picasso Disposable Laser Tips

Variety Kit (25 Assorted Tips): TIP25

20 Endo Tips (200µm 20mm): TIP2006020

25 Surgical EZ Tips (300µm 5mm): TIP3006005

25 Surgical Tips (300µm 10mm): TIP3006010

25 Standard Tips (400µm 5mm): TIP4006005

25 Perio Tips (400µm 10mm): TIP4006010

Stripped, cleaved, and ready-to-go, disposable laser tips are the easiest and fastest way to transition between patients and dental procedures.



Replace or resupply your laser tips with our patented ready-to-use disposable or reusable tips for your diode laser.

Strippable Fiber

SF-KIT-PLUS

All your Picasso laser strippable fiber tools, one convenient kit. Each kit contains 1 fiber stripper, 1 cleaving stone, 1 spool, 2 surgical handpieces, and 20 disposable tips. Available in PLUS or FC connections.



Picasso Fiber Spool

MTHPFS

Fiber Spool for Multi-Tip Handpiece to use with all models of Picasso lasers that use disposable tips, NOT with strippable fiber.



Multi-tip handpeice and fiber cable replacement for Picasso Series Lasers.



002-00390

60° Bend disposable tips that are made specifically to fit your surgical handpiece. These tips snap right on to your surgical handpiece for a perfect bend of your strippable fiber.



Picasso Whitening Quadra Tip

QT

This Quarda Tip is used laser assisted whitening with the Picasso Plus. It is compatible with the multi-tip hand piece.

THE FASTEST AND MOST CONSISTENT CURE IN DENTISTRY



YOUR PATIENTS ARE READY FOR LIGHTSPEED.


Patients are ready for improved care through the latest in dental technology; they are actively searching for practices with the most up-to-date and innovative patient care experiences. Soft cures can leave patients with pain, prolonged sensitivity, and a weak bond can open the tooth to future structural damage. Nothing disrupts a practice like an incomplete cure, with proper use of the Monet, this is now a problem of the past.

Monet Laser Curing Light 001-00105

FROM THE INNOVATIVE MIND THAT INVENTED LED CURING COMES A QUANTUM LEAP IN CURING TECHNOLOGY: THE WORLD'S FIRST HANDHELD LASER CURING LIGHT.

"High irradiance output produced faster and deeper polymerization than other light sources"

Clinicians Report | Gordon J. Christensen

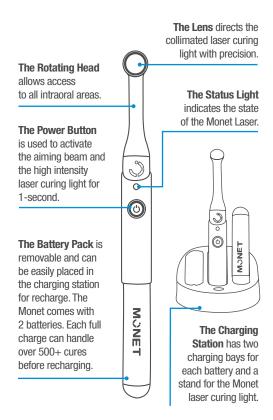

The handheld Monet curing laser is the first of its kind. The collimated beam and consistent power of the Monet laser makes for superior bond strength, and faster, deeper, more reliable curing results.

LASER CURING **MEANS CONFIDENCE.**

Laser curing has major advantages over LED curing. With proper application of a Monet laser cure, you will never have to worry about incomplete bonds or soft cures again. The Monet laser produces a deeper cure than LED or any other light sources during the same exposure time.

LASER CURING **MEANS CONSISTENCE.**

Laser curing provides a consistent dispersion of energy and intensity at any distance to create a complete, reliable, and even beam pattern. The collimated beam improves polymerization through composites. Each click is a one-second cure reaching through 2.5mm of composite depth, that means it's possible that bulk fills can be done in three seconds!



Valo is a registered trademark of Ultradent Products; Sapphire is a registered trademark of DenMat, Kerr is Trademarked to Demetron®

DEPTH OF CURE










It is generally recommended to do a single 1-second cure with a composite layer less than 2.5 mm, two 1-second cures for layers between 2.5-5 mm, and three 1-second cures for layers more than 5 mm. To minimize heat, do not cure any spot more than 3 seconds in rapid succession. For prep areas larger than 8 mm in diameter, two overlapping spot cures is recommended. Monet comes with unique apertures for precision control of your laser. For preps near gingival tissues, consider using small aperture of 2, 4, or 6 mm to avoid light on gingival tissue. It also comes with a reducing aperture which cuts the power by 50% to give you even more custom control!



World's First Take-Home and In-Office PROFESSIONAL TEETH WHITENING STRIP SYSTEM Available ONLY From Your Dentist

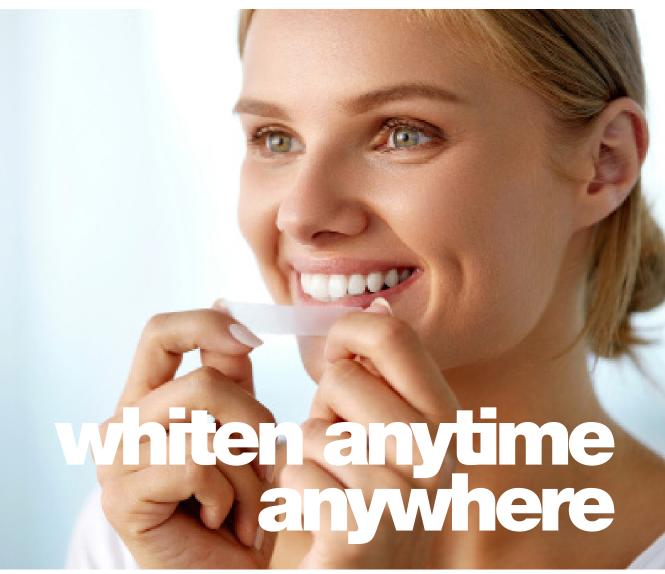

Sheer White!® professional whitening strips stick and mold perfectly to the teeth, becoming almost invisible. They're comfortable to wear during any activity and cause little-to-no-sensitivity.

Exercise

Offer your patients a powerful teeth whitening solution that is fast, effective, and virtually painless. Sheer White! allows you to offer whiteing to all your patients by giving you a wallet friendly option to other professional whitening products. Custom trays and gels are often too expensive for all of your patients and can be painful to some. Sheer White! will save you money and create repeat whitening visits. Your patients will love their whiter smile as much as they love saving time and money. **Better Revenue, Better Whitening!**

In-Office



Step 1 Sheer White! In-Office

Apply an In-Office whitening strip at the end of any appointment, application onlyn take 1-minute. The SheerWhite! In-Office strip is a 20% Hydrogen Peroxide formula delivered by our patented non-slip strip, and works with minimal sensitivity. SheerWhite In-Office is capable of whitening up to 3 shades whiter in just 30-60 minutes.

PATIENT PACK OF 2 006-00234 6 PATIENT KIT 006-00251

Take-Home



Step 2 Sheer White! Take-Home

Send your patients home with the Sheer White!
Take-Home patient treatment kit, a 20% carbamide peroxide formula delivered by our patented go-anywhere non-slip whiteing strip.
Each SheerWhite! Take-Home kit includes a 5-day treatment, consisting of 10 strips. Patients will see up to 5 more shadesof whitening in just 5 days.

6 PATIENT INTRO KIT 006-00081 24 PATIENT BULK PACK 006-00099

Sheer White In-Office and Take Home can be used on thier own, but for best results and patients who will see up to 8 shades whiter, use as a system. Apply the In-Office strips at the end of any appointment then send your patients on their way with a **whiten anytime, anywhere** Sheer White! take-home kit.



12 BOTTLES 006-00266 48 BOTTLES 006-00267

Cavity pHighter Mints

12 BOTTLES 006-00264 48 BOTTLES 006-00265

Two critical factors contributing to the risk of dental caries are pH level and buffer capacity in saliva. Cavity pHighter strips can quickly and accurately identify both to rapidly indicate a risk of caries. It's so easy to use that you can test in office or send the strips home with a patient to track their risk.

Improving a patients hygiene habits is easy with Cavity pHighter. The mints increase saliva pH and saliva flow, as well as bolster buffer capacity by releasing soluble calcium all to promote healthy tooth enamel.

- Create additional revenue from salivary testing.
- Increase revenue by supplementing individual care plans
- Continued revenue from providing test strips and mints for home use.
- Generate patient engagement and retention



AO GROUP INC



Lasers& Curing& Whitening& Smiles.

